<i>Kojenecká váha s aplikací • NÁVOD K OBSLUZE</i>	<b>CZ</b>	2-7
<i>Dojčenská váha s aplikáciou • NÁVOD NA OBSLUHU</i>	<b>SK</b>	8-13
<i>Baby scale with app • USER MANUAL</i>	<b>EN</b>	14-19
<i>Csecsemő mérleg alkalmazással • HASZNÁLATI ÚTMUTATÓ</i>	<b>HU</b>	20-25
<i>Waga niemowlęca z aplikacją • INSTRUKCJA OBSŁUGI</i>	<b>PL</b>	26-32
<i>Babywaage mit App • BEDIENUNGSANLEITUNG</i>	<b>DE</b>	33-39

# Mimi

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**INSTRUCTIONS FOR USE**

Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the receipt and, if possible, the box with the internal packing.

**Tip**

- For initial pairing with the "AiLink" app and instructions on how to use the app, please visit our website [www.etasince1943.com/baby](http://www.etasince1943.com/baby).

**I. SAFETY WARNING**

- Consider the instructions for use as a part of the appliance and pass them on to any other user of the appliance.
- The appliance is intended for use in households only and for similar use! It is not intended for use in health-care facilities and for commercial use!
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Never use the appliance if it malfunctions, if it fell to the ground and got damaged or if it fell into water. In these cases bring the appliance to a specialised service to have the safety and functions of the device checked.
- **Never immerse the scale into water (even their parts)!**
- If the scale was stored in lower temperatures, let it adjust to the room temperature.
- Use the scale in a horizontal position only and at places where they cannot knock over, in sufficient distance from heat sources (e.g. **stove, fireplace, heat radiators**), sources of increased air humidity (e.g. **sauna, bathroom, swimming pool**) and equipment with strong electromagnetic field (e.g. **microwave oven, radio**).
- Do not place the scale on uneven surface or on surface covered with carpet. Unstable or soft surface under the scale can have negative effect on precision of weighing.
- In handling the scale, make sure injury of persons or damage of the appliance is prevented.
- Do not overload the appliance weight greater than the maximum weight capacity! This prevents damage to the appliance.
- Do not put any objects on the surface of the scale during storing, they could be damaged.
- Do not use rechargeable accumulators.
- If the battery is leaking, remove it immediately; it can damage the scale.
- Dispose of a discharged batteries in a suitable way (see par. **VIII. ENVIRONMENTAL PROTECTION**).

- Keep batteries and scale out of reach of children and incapacitated people. A person, who would swallow the batteries, must seek medical help immediately.
- Never use the appliance for any other purpose than for that described in these instructions for use!
- **WARNING:** There is a risk of injury in case of incorrect use (not in accordance with the manual) of the appliance.
- The producer is not responsible for damage caused by improper use of the appliance (e.g. **breaking the board**, etc.).

#### Important notes:

- Exercise extra caution when handling the baby, always hold the baby and keep a close watch.
- We also recommend placing a soft pad (**towel, diaper**, etc.) on the weighing surface to keep the baby comfortable.

#### Wireless connectivity and compatibility

- ETA can in no way be held liable for any accidental, indirect or subsequent damage, data loss or any other loss incurred as a result of information leakage during communication through BLUETOOTH, Wifi technology.
- 100% successful wireless communication cannot be guaranteed with all Smart devices. Devices have to be compatible with BLUETOOTH specifications set forth by the company Bluetooth SIG, Inc., and their compliance with the specifications has to be certified. However, given the variety of devices available in the market, even if a device is in compliance with the specifications, there might be cases when the character or specifications of the given Smart device do not allow connectivity, or hinder connectivity in any way. This is not a reason to make a warranty claim.
- Radio waves transmitted by devices may interfere with operation of medical equipment. Such interference may cause a failure. Therefore, do not use the device in the following places:
  - in hospitals, trains, planes, petrol stations and in areas where flammable gases may be present.
  - in the vicinity of automatic doors or fire detectors.

**The application is owned and operated by a third party (hereinafter the “Subject”) different from HP TRONIC Zlín, Ltd. (hereinafter the „Seller“). No personal information is transmitted between the Subject and the Seller. When registering for the application, the customer transfers the personal data directly to the Subject, which becomes the personal data controller. The Application and the associated privacy policy are governed by the Rules of the Subject. Development and any changes to the application are made solely by the owner of the application.**

#### Note

- The application requires an active internet connection.

## II. GENERAL INFORMATION

After the scale is paired with the Smart device, the AiLink app starts displaying the measurement results in real time. The scale itself shows only the weight, but the application allows you to enter other values, with which it can further work. Before first use, download the AiLink application to your Smart device for free. You can find the application in the Apple Store or Google Play database, or to speed up the process, you can use the following QR codes:

### Note

- If the QR code is not active, find the application manually by name.

Detailed instructions and info on the AiLink app itself can be found on the website:

[www.etasince1943.com/baby](http://www.etasince1943.com/baby)

## III. PREPARATION FOR USE

Remove all the packing material and take the scale out. Remove all possible adhesion foils, stick-on labels or paper from the scale. Open the battery cover on the bottom of the scale, **insert the battery with the correct polarity** and close the cover in the opposite direction. Place the scale on flat, stable and dry surface.

### Units

The scale itself has a weight unit switch (the **A1 "UNIT"** button), so you can change the units between g - kg - lb - st: lb. Weight units can also be changed in the application.

kg = kilogram	1 kg = 2.2 lb
lb = pound (British weight unit)	1 lb = 0.454 kg
st = stone (American weight unit)	1 st = 6,356 kg

### Replacing batteries

Open the battery cover on the bottom of the scale. Insert a battery with the right polarity and close the cover. Use AAA batteries (3 pcs) for this appliance. If you are not using the appliance for some time, take the batteries out!

### Notes

- Do not combine various types of batteries. Do not even combine used and new batteries together.
- Batteries not intended for this purpose must not be charged.
- Remove the batteries if you will not be using the appliance for a long time.
- The supply terminals must not be short-circuited.

## IV. INSTRUCTIONS FOR USE

For the most comparable measurement results, ideally weigh the baby under the same conditions and at a similar time. Feeding, clothing, etc. can make a big difference between measurements in small babies. Short-term weight fluctuations are generally attributed to fluid depletion (loss).

**UNIT button (A1)** - used for switching weighing units.

**ON/OFF // TARE button (A2)** - press it in the off state to switch the scale on. Press during weighing to activate TARE mode and long press to deactivate the scale.

### 1) STANDARD WEIGHING

Without using the AiLink app on the Smart device, the device provides the basic weight weighing function found on all conventional scales.

- Switch the scale on by pressing the “A2” button. The display will show “HI”.
- It is possible to place, for example, a blanket on the weighing surface so that the baby lies on a soft surface. Then press the “A2” button to reset the weight of the blanket (TARE function).
- When the display shows 0, carefully place the baby on the weighing surface.
- Try to keep the baby as calm as possible. After a moment, the weighed weight will flash 3 times and the weight value will lock. Then the scale will automatically switch off after a while.

### Notes

- The scale is switched on automatically if you do not step on it within approx. 30 seconds after „0.0“ is displayed.
- The icon flashes on the display to indicate communication / searching for nearby devices. When the icon stays lit, the connection with the application is on.
- If the weight is more than 3 kg and you press the on/off button for the TARE function, the display will reset and “TARE” will be displayed at the same time.

## 2) WEIGHING USING THE APP

Download the app, create an account and pair the app with the scale. Turn the app on and then weigh the baby in the standard way. Once successfully weighed, the data will automatically be transferred to the app, where you can add additional data to the weight and then track it within the weighing history.

## V. TROUBLESHOOTING OF POSSIBLE PROBLEMS

### If you have problems with weight:

- Check whether the battery is inserted properly.
- Check whether the scale is on a flat pad and they do not touch the wall.
- If the scale cannot be switched on or only “LO” is displayed, replace the batteries.
- If „Err“ appears on the display, the scale was overloaded.
- If unusual signs are displayed, remove the insert the battery again.
- If signs appear, such as if the display is flashing or if there is a faulty sign, try to move the scale away from the influence of a source of interference or ensure switching off the source for the time of using the scale.

## VI. STORING

After using the scale, put it into the horizontal position (not vertically), out of reach of children and incapacitated people.


## VII. MAINTENANCE

Do not use coarse and aggressive cleaning agents (e.g. **sharp objects, scrapers, diluting agents or other solvents**). Clean the surface of the scale surface by wiping it with a soft damp cloth. Make sure that water does not get inside the appliance.

**More extensive maintenance or maintenance requiring intervention in the inner parts of the appliance must be carried out by professional service!**

**Failure to comply with the instructions of the manufacturer will lead to expiration of the right to repair!**

The information regarding warranty and product repairs you can find on the website [www.eta.cz/servis-eu](http://www.eta.cz/servis-eu).
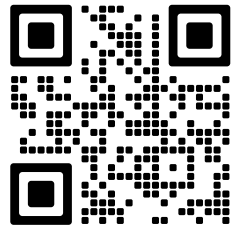
## VIII. ENVIRONMENTAL PROTECTION

If the dimensions allow, there are labels of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal, hand them over at special collection places where they will be accepted free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal.

Ask for more details at the local authorities or at a collection site. Fines can be imposed for improper disposal of this type of waste in agreement with the national regulations. Remove a discharged battery and dispose of it in a suitable and safe way through the special collection networks. Never dispose of the battery by burning!

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Weighing range max.:	20 kg / 44 lb / 3,1 st
Discrimination:	5 g
Power supply:	3 x Batteries AAA (included)
Net weight of device:	2,1 kg
Size of the product (mm):	630 x 300 x 54
	Protection class of the appliance: III.
Frequency range:	2402 - 2480 MHz
Output power (max):	0,3 dBm
Bluetooth version:	4.2
Compatible operating systems:	Android 6.0; iOS 9.0 (and higher)
Maximum communication reach:	10 m
	(open space without any obstacles and interferences)

Hereby, ETA a.s. declares that the radio equipment type ETA7306 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: <http://eta.cz/declaration-of-conformity>

**The manufacturer reserves the right to change the technical specifications for the respective model. The manufacturer reserves the right to change/update the AiLink application without prior notice. Due to the character of such changes/ updates, some information or procedures in this manual might become obsolete or might be missing.**

### WARNINGS AND SYMBOLS USED ON THE APPLIANCE, PACKAGING OR IN THE INSTRUCTIONS MANUAL:

HOUSEHOLD USE ONLY. DO NOT IMMERSE IN WATER OR OTHER LIQUIDS.	
	TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, ARRANGES OR PLAYPENS. THIS BAG IS NOT A TOY.
	CAUTION
	Please read the instructions manual

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