

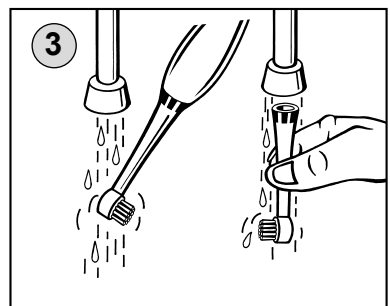
<i>Bateriový zubní kartáček • NÁVOD K OBSLUZE</i>	<b>CZ</b> 4-7
<i>Batériová zubná kefka • NÁVOD NA OBSLUHU</i>	<b>SK</b> 8-11
<i>Battery toothbrush • INSTRUCTIONS FOR USE</i>	<b>GB</b> 12-15
<i>Elektromos fogkefe • HASZNÁLATI UTASÍTÁS</i>	<b>H</b> 16-19
<i>Szczoteczka do zębów na baterie • INSTRUKCJA OBSŁUGI</i>	<b>PL</b> 20-23

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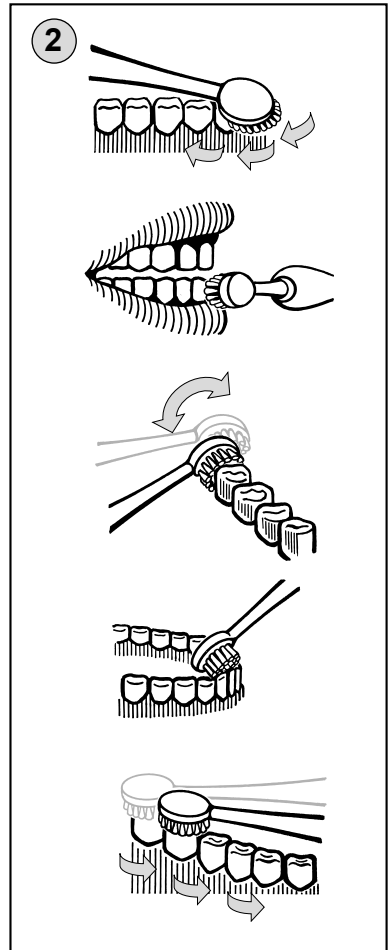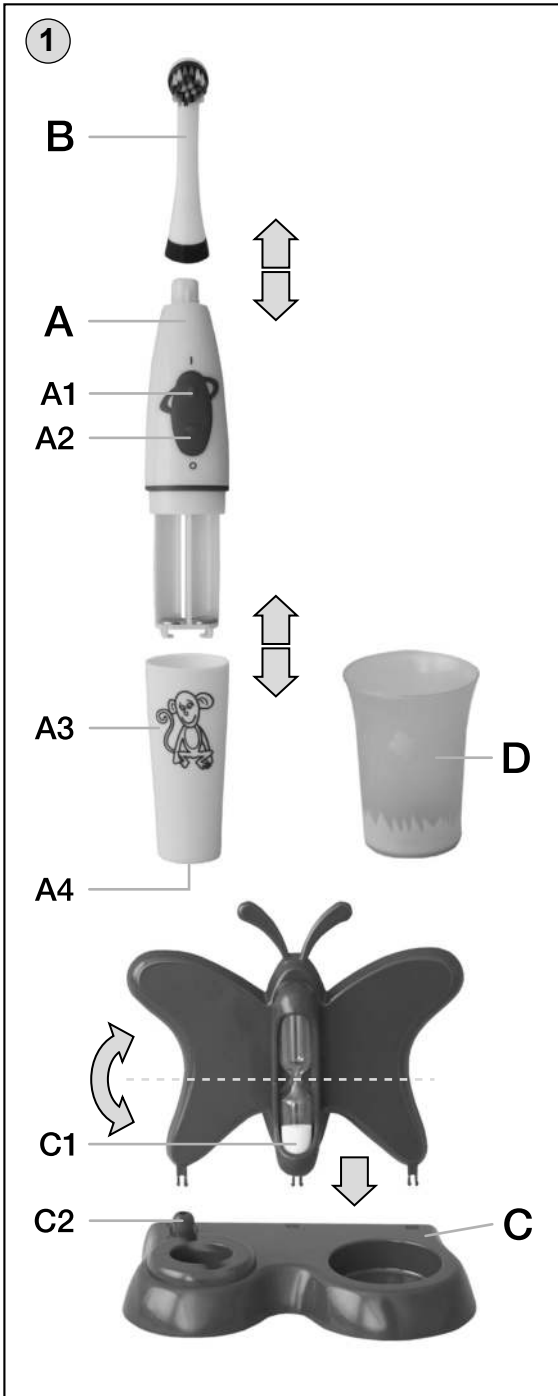# eta 1294

## USER MANUAL

Dear customer, thank you for purchasing our product. Prior to putting this device into operation, please read these instructions for use carefully and be sure to store well these instructions along with the warranty card, receipt and, if possible, also the product packaging and its inner contents. This device is a mechanical remover of tooth plaque which is often used and which will help your children improve their oral hygiene. Tooth brush represents the latest technology contained in a effective and ergonomic device which is an ideal tool for teeth care.

## I. GENERAL SAFETY INSTRUCTIONS

- Consider the instructions for use as a part of the appliance and pass them on to any other user of the appliance.
- Never use the device if it's not working properly, has fallen on the ground or has been damaged.
- This product is not intended for children under 3. Children at the age of 3 to 14 and persons with limited physical, sensory or mental capabilities may only use the toothbrush if they are supervised or have been instructed of the use of the appliance in a safe manner and understand possible risks. Children must not play with the appliance. Children must not perform cleaning and maintenance intended to be performed by users, unless they are older than 8 years and are under supervision. Children under 8 years of age must be kept of out reach of the appliance.
- The appliance is only intended for individual use at home for personal oral cavity care! It is not intended for commercial use, for more than one patient
- Never immerse the propulsion unit in water and other liquids!
- Use the appliance only with accessory (replacement brush heads) intended for this type of device.
- Do not put the brush away on hot heat sources (e.g. **stove, radiator, cooker**)
- If the device has been stored at lower temperatures, acclimatize it first which will maintain the mechanical compactness of the plastic parts.
- Do not insert the appliance or its accessories in any body cavities, except for oral cavity.
- The appliance is not designed to care for animal teeth.
- Clean the accessories after each use. Failure to adhere to this instruction might cause worsening your state of health, because without disinfection, spread of mildew and bacteria is supported
- Do not leave the appliance unattended!
- Protect the appliance from dust, chemicals, toiletries, liquid cosmetics, large fluctuations in temperature, and direct sunshine
- Clean the brush after use (see par. **IV. MAINTANENCE**) and store it at a dry place.
- Do not use the appliance for any other use than the one it is intended for and described!
- The manufacturer is not liable for any damage caused by incorrect use of the device (e.g. injury of oral cavity, gums, teeth, damaging the device) and is not liable for the warranties in case of breaching the precautions mentioned above.

## SPECIAL SAFETY INSTRUCTIONS

- If you have unpleasant feeling in your mouth or if you feel pain during the cleaning, stop using the appliance and contact your dentist. Contact the dentist also in the case that you have any doubt about use of the appliance in your individual case.
- If your toothpaste contains peroxide, baking soda or other bicarbonate (frequently used in whitening toothpastes), thoroughly clean the brush and its handle after use to prevent damaging the appliance.
- Do not use the appliance in care for incapacitated persons or persons who are unable to express their feelings.
- **If you have undergone tooth or gum operation recently, suffered long-term problems with your teeth or were subject to extensive medical interventions in the last 2 months, or if you suspect problems of dental character (with fillings, etc.), consult your dentist before using the appliance.**
- **If you experience any health troubles (for example considerable bleeding, or if any bleeding persists even after 2 weeks of use of the brush), stop using the appliance and contact your physician for consultation.**
- For hygienic reasons do not share the toothbrush with another person or more persons. Infections could spread or oral cavity inflammation could be initiated.
- Only use the appliance with the original accessories designated for the respective type. Do not use damaged/deformed adapters.
- **WARNING:** There is a risk of injury in case of incorrect use (not in accordance with the manual) of the appliance.
- All text in other languages and pictures shown on the packaging, or the product itself, are translated and explained at the end of this language mutation.

## SPECIAL SAFETY INSTRUCTIONS

- If you are not going to use the appliance for a long time, remove the batteries.
- After the batteries are completely empty, replace them immediately with new ones.
- In the case of battery leak, immediately replace the battery, otherwise it may damage the appliance.
- Do not mix the old and new batteries.
- Do not use rechargeable batteries.
- Keep the batteries out of reach of children and incapacitated persons. If you swallow the battery, immediately seek for medical assistance.
- Dispose the used batteries in a proper manner (see V. ENVIRONMENTAL PROTECTION section).

## II. FEATURES OF THE DEVICE AND ACCESSORIES (pic. 1)

### A – propulsion unit

- A1 – button I (switch)
- A2 – button 0 (switch)

- A3 – battery cover
- A4 – battery cover arrest control

### B – brush head

### C – stand

- C1 – timer

- C2 – brush head holder

### D – cup

## III. INSTRUCTIONS FOR USE

Remove all packaging material and remove the toothbrush with accessories. Rinse the parts that will come in touch with the oral cavity in clean hot water before first use.

### Notes

- You may use any quality tooth paste for cleaning your teeth. We do not, however, recommend using special tooth pastes with whitening effects. In order to prevent spraying, first put the brush on the teeth and then turn on the propulsion unit.
- Do not apply excessive pressure on your teeth.
- You can also use this brush when you have fixed braces, facets, crowns, etc. However, please bear in mind that if the brush is used with braces, the brush bristles are subject to more wear.
- We recommend that you replace the brush (**B**) with a new one after 3 months of use. However, if the brush is showing signs of significant wear or damage.

### Inserting batteries

Release the battery cover **A3** on the bottom of the propulsion unit **A** by turning the control **A4** clock-wise into a position ● and remove the cover. Insert 2 pcs. of AAA batteries, observing the +/- signs (see symbol). Replace the cover and arrest it by turning the control into the position ■.

### Use

Slide the brush head **B** slightly on the propulsion unit **A** so that its protrusions fit onto the protrusions of the propulsion unit. Rinse your mouth thoroughly before cleaning your teeth, put toothpaste on the brush head and moisten it. Position the brush against the teeth so that it covers the surface of the teeth as well as the gum line. Turn the brush on by pressing the **I** button. Guide the brush slowly along the teeth and grinders, follow the curve of the gums and the shape of the individual teeth. Hold the brush always on one spot for several seconds before you move over to a next tooth. Make sure you have cleaned the whole surface of the teeth (the inner and outer part of the occlusion area), as well as the line around the gums (pic. 2). Do not press on the brush, leave it to do all the work by itself. Acquired and permanently maintained procedure will prevent you from omitting some parts of the teeth. Do the cleaning for at least 2 minutes which is the minimal recommended time by stomathologists. You may use the timer **C1** for this purpose. If the plaque is removed properly, the gums will gain pink colour and will surround the teeth fast.

### Timer (C1)

The timer in the form of hourglass allows to control the length of teeth cleaning (approx. 2 min.).

## IV. CLEANING AND MAINTENANCE

**Clean the toothbrush regularly after each use! Do not use any rough and abrasive cleaning agents!** Wipe the surface with a soft damp cloth. Rinse the brush head thoroughly for a few seconds under running water while leaving the propulsion unit on (pic. 3). Turn the propulsion unit by pressing the **0** button, take off the brush head and rinse both the parts in running water separately under running water (pic. 3) and let dry. To achieve hygiene cleanliness we recommend immersing the brush head every 2 - 4 in a cleaning agent for dentures.

### Replacing the brush head

To achieve the best results while cleaning teeth, it is necessary to replace the brush head every time its bristles start to bend outward, i.e. after approx. 3 months. You can buy the replacement head at your distributor.

## V. ENVIRONMENTAL PROTECTION

If dimensions permit, all pieces have printed symbols of materials used for production of packing, components and accessories as well as recycling symbols. The symbols on the product or documentation means the electric or electronic parts used in the appliance may not be disposed off together with communal waste. Hand over the electric or electronic part to a waste collection yard for free disposal. Proper disposal of the product will help preserving valuable natural resources and preventing potential negative impacts on environment and human health that would otherwise result from improper waste disposal. For more details, please contact your local authorities or collection yard. Fines may be imposed for improper disposal of this type of waste pursuant to national directives. Remove flat batteries from the appliance and dispose them off in a suitable way using special collection centres. Do not incinerate the batteries!

**Specialized service is required for overhaul maintenance or maintenance requiring disassembly of the internal components of the appliance! The right for warranty repair expires if the above mentioned instructions are not adhered to.**

## VI. TECHNICAL DATA

Weight (kg) about 0,1  
 Noise level: Acoustic noise level of 55 dB (A) re 1pW

**The manufacturer reserves its right for unimportant deviations from the standard design without influence of the product function.**

**WARNINGS AND SYMBOLS USED ON THE APPLIANCE, PACKAGING OR IN THE INSTRUCTIONS MANUAL:**

*HOUSEHOLD USE ONLY.  
 DO NOT IMMERSE IN WATER.*

*TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.*

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© DATE 16/7/2013