beurer

BF 100 Body Complete

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ENGLISH

 Comprehensive instructions for operating the device are provided in the "Instructions for use" PDF file from the download area at <http://www.beurer.com>.

1. Safety notes

Safety notes

- The scales must not be used by persons with medical implants (e.g. heart pacemakers). Otherwise their function could be impaired.
- Do not use during pregnancy.

- Do not stand on the outer edge to one side of the scale otherwise it may tip!
- Caution! Do not step onto the scale with wet feet or when the scale's surface is damp – danger of slipping!
- Keep children away from packaging materials (risk of suffocation).

General notes

- The unit is for personal use only and is not intended for medical or commercial applications.
- Place the scale on a firm level floor; a firm floor covering is required for correct measurement.

- Repairs may only be carried out by customer services or authorised retailers. Before submitting any complaint, first check the batteries and replace them if necessary.
- Should you have any questions concerning the use of our devices, please contact your dealer or customer service.

Storage and maintenance


The accuracy of the measured values and the service life of the unit depend on careful handling:

Important

- The unit should be cleaned occasionally. Do not use abrasive detergents and never immerse the unit in water.
- Ensure that no liquid gets on the scale. Never dip the scale into water. Never wash it under running water.
- Protect the scale from impacts, moisture, dust, chemicals, large changes in temperature and sources of heat that are too close (oven, radiator). Do not press the buttons using excessive force or with pointed objects.

2. Description of the unit

Overview

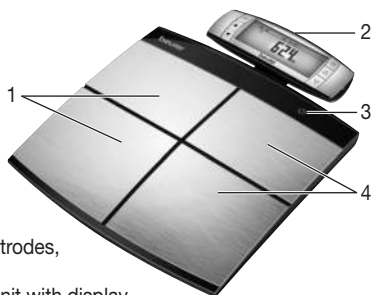
- 1 Foot electrodes, left foot
- 2 Control unit with display
- 3 Operating status indicator
- 4 Foot electrodes, right foot

Control unit with display

-
- Diagram of the control unit with display showing buttons and indicators labeled 5 through 12.
- 5 Contact electrodes, left hand
 - 6 Contact electrodes, right hand
 - 7 set Setting button, confirms inputs
 - 8 Menu: up button
 - 9 Menu: down button
 - 10 ON button
 - 11 Settings: up button
 - 12 Settings: down button

The display

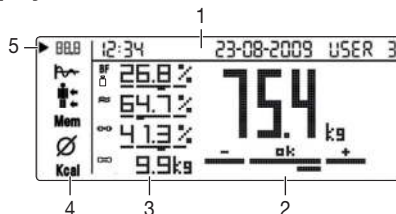
- 1 Time, date, user
- 2 Weight with analysis
 Weight
- 3 Analysis data (partly with assessment)
BF Body fat
 Water
 Muscle percentage
 Bone mass
- 4 Menu
 Overview of current measurement
 Trend
 Detailed display of body fat, muscle percentage of upper/lower body
Mem Saved measurements 1–30
 Average values
Kcal Calorie consumption
- 5 Indicator for the current menu item


3. Initial use

Inserting batteries

Control unit

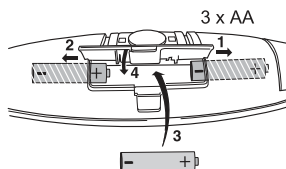Scale

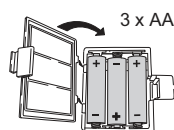Basic settings

a. Setting the unit

[SET] → "SET_UNIT"

Sequence	Set values
Language	
Measurement unit	kg – lb – st
Time	24 h – 12 h
Guest mode	On – Off

b. Setting the time

[SET] → "SET_TIME"

Sequence	Set values
Date	Year – Month – Day
Time	Minutes – Hours

Use the [] and [] buttons to select and confirm with [SET].

Setting user data

[SET] → "SET_USER"

User data	Set values
Sex	male  , female 
Height	100 to 250 cm (3' 03" to 8' 02")
Date of birth	year, month, day
Degree of activity	1 to 5
Target weight	adjustable, in kg (for graphical progress display)







Use the [▽] and [△] buttons to select and confirm with [SET].

Degrees of activity

Degree of activity	Physical activity
1	None.
2	Low: A small amount of light physical effort (e.g. short walks, light garden work, gymnastic exercises).
3	Medium: Physical effort for 30 minutes at least 2 to 4 times a week.
4	High: Physical effort for 30 minutes at least 4 to 6 times a week.
5	Very high: Intensive physical effort, intensive training or hard physical work for at least one hour daily.


4. Weighing with diagnostics

[P] → 0.0 kg → [▽]/[△] [SET]

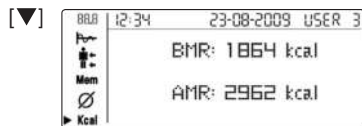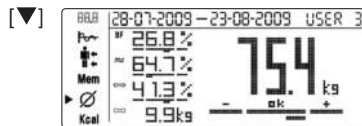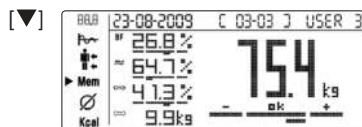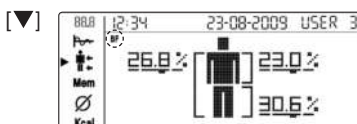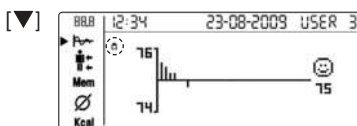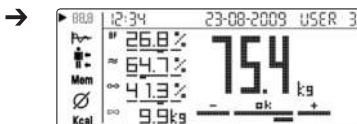5. Evaluation of results

The measurements are shown on the display and evaluated using the bar chart for simple classification:

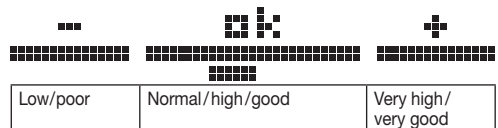6. PC transfer

System requirements:

Windows 7, Windows Vista, Windows XP SP3

HealthManager:


Health management in conjunction with Beurer blood pressure and blood glucose monitors.

EasyFit:

Training analysis in conjunction with Beurer heart rate monitors.

Process:

- You will find the Beurer software download and installation in the download area at <http://www.beurer.com>
- Establish USB connection.

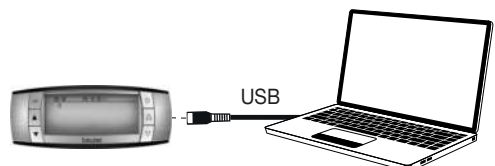- Follow the instructions for software installation.

7. Disposal

Standard and rechargeable batteries should not be disposed of separately from the household waste. As a consumer, you are legally obliged to return used batteries for proper disposal. You can hand in your used batteries at

public collection points in your district or sales outlets where batteries of this type are sold.

Note

The codes below are printed on batteries containing harmful substances:

Pb = battery contains lead,

Cd = battery contains cadmium,

Hg = battery contains mercury.

For environmental reasons, do not dispose of the scale with the batteries in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point.

Observe the local regulations for material disposal.

Dispose of the unit in accordance with EC Directive 2002/96/EC-WEEE (Waste Electrical and Electronic Equipment).

If you have any questions, please contact the local authorities responsible for waste disposal.

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